

Hydrophilic Nelaton with Water Sachet Coude Tip

DESCRIPTION

m|devices Hydrophilic Nelaton Catheter with Water Sachet Coude Tip contains a blue coloured holding sleeve for easy visualisation for an aseptic insertion. This all in one product offers convenience and ease of use without the need for additional products. The hydrophilic coating in conjunction with the Coude angled closed distal end and bilateral side eyelets, support a smooth and painless insertion and efficient drainage of urine.

FEATURES & BENEFITS

- Hydrophilic coated
- All-in-one product with water sachet
- Blue coloured holding sleeve to maintain sterility during insertion
- Coude tip to support difficult insertions
- Sterile, easy peel packaging with adhesive wall attachment
- Colour coded connectors
- Heat polished eyelets
- Smooth round tip


PRODUCT DATA

- Length: 40cm
- Sizes 8Fr-18Fr
- Not Made with Natural Rubber Latex
- Not Made with DEHP
- Material: Nelaton: PVC
- Sachet: Purified water
- Holding sleeve: PE
- Shelf life: 3 years
- ARTG number: 210992

ISO STANDARDS

- EN ISO 14155-1:2003 Clinical investigation for medical devices for human subjects – part1 General requirements
- EN550:1994 Sterilisation of medical devices – Validation and routine control of ethylene oxide sterilisation
- EN ISO11607-1:2006 Packaging for terminally sterilised medical devices – Part 1: Requirements for materials, sterile barrier systems and packaging systems

PRODUCT DIAGRAM

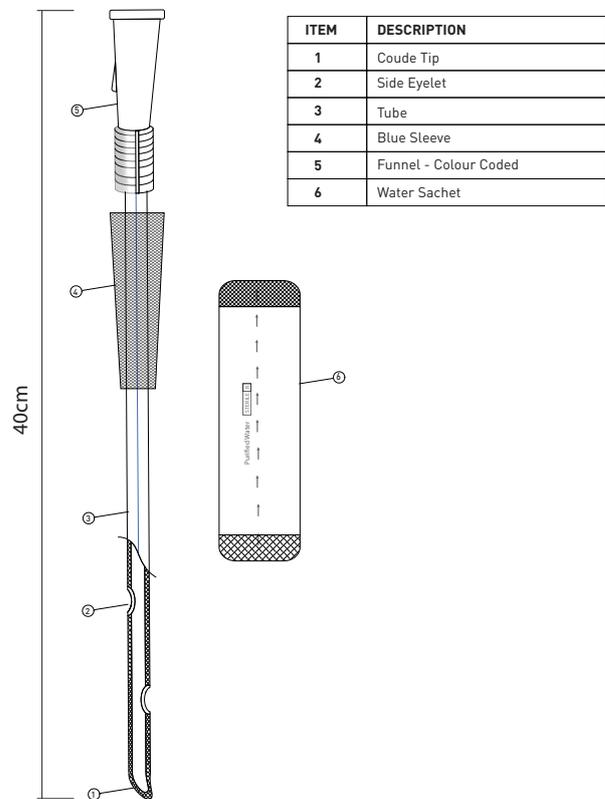

ITEM	DESCRIPTION
1	Coude Tip
2	Side Eyelet
3	Tube
4	Blue Sleeve
5	Funnel - Colour Coded
6	Water Sachet

Hydrophilic Nelaton with Water Sachet Coude Tip

Hydrophilic Nelaton with Water Sachet					STERILE
RE-ORDER	DESCRIPTION	LENGTH	GAUGE	COLOUR	CARTON/QTY
UR049001	Hydrophilic Nelaton with Water Sachet - Coude Tip	40cm	10Fr	Black	25
UR049002	Hydrophilic Nelaton with Water Sachet - Coude Tip	40cm	12Fr	White	25
UR049003	Hydrophilic Nelaton with Water Sachet - Coude Tip	40cm	14Fr	Green	25
UR049004	Hydrophilic Nelaton with Water Sachet - Coude Tip	40cm	16Fr	Orange	25
UR049005	Hydrophilic Nelaton with Water Sachet - Coude Tip	40cm	18Fr	Red	25